

All-body Cushions

Functions:

1、All-In-One board (including extension component)

The all-in-one board can be compatible with all cushions, and also can be used as supine head-neck-shoulders baseboard and belly-pelvis baseboard.

2、 Prone Pelvic Cushions

It has three sizes of padded with hole, there are small one, medium one and large one, and one padded without hole inserts to accommodate a variety of patients. Load the cushions on the base plate directly and used for prone pelvis fixation and craniospinal irradiation fixation.

3、 Prone Breast Cushions


The cavity wedge for the breast is available for the contralateral breast to assist in proper positioning. The cushions are used for breast conservation therapy in prone position and fixation.

4、 Supine Breast Cushions

There are three flat cushions, which heights are 70mm, 85mm and 101mm; three wedge-shape cushions, which angles are 5°, 10° and 15°. The cushions are used for breast therapy in supine position and fixation.

5、 Lower Limbs Cushions

The lower limbs positioning system include one carbon fiber base plate, two knees cushions and two feet cushions and used for the lower limbs therapy. This system can be used as lower limbs fixation together with all-in-one board and all other cushions, such as prone pelvis cushions, prone breast cushions, supine breast cushions, and the knees cushions used as ankle fixation in prone position.

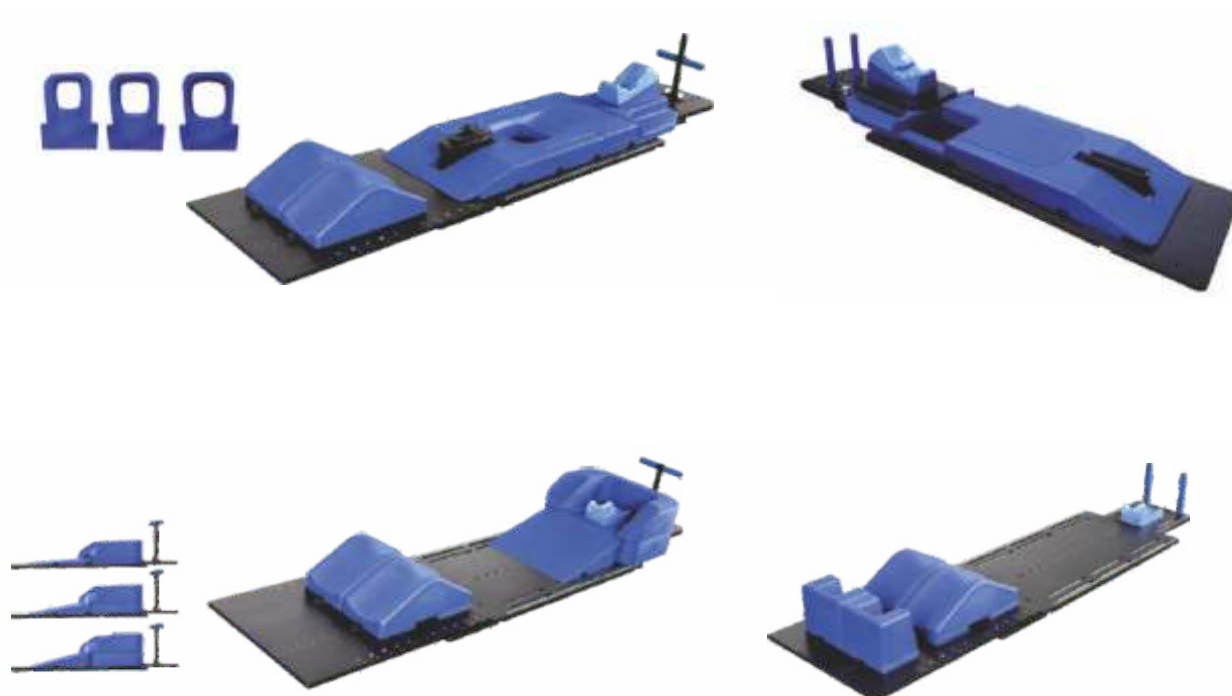

Supine Breast Cushions



YC-KIV-155

1. Polyurethane material;
2. Load it on the base plate directly, firmly fixed and repetitive used;
3. Supported unilateral breast treatment;
4. Three angles of three wedge-shape cushions: 5°, 10°, 15°; and the heights of three plat cushions: 70mm, 85mm and 101mm;
5. Ergonomic hand grips and arm cushions comfortably support the patient.
6. Not including the All-In-One board.

Lower Limbs Cushions




YC-KIV-158

1. Polyurethane material;
2. Can be used for single lower limbs therapy fixation;
3. Can be used as lower limbs fixation in craniospinal irradiation therapy;
4. Not including the All-In-One board.

All-In-One Board


All-In-One Board (T-type) >>>

YC-KIII-136T

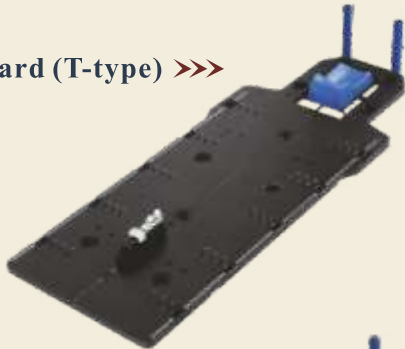

All-In-One Board >>>

YC-KIII-136P

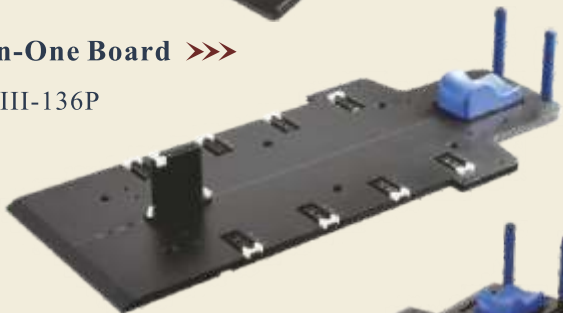

All-In-One Board >>>

YC-KIII-135

